
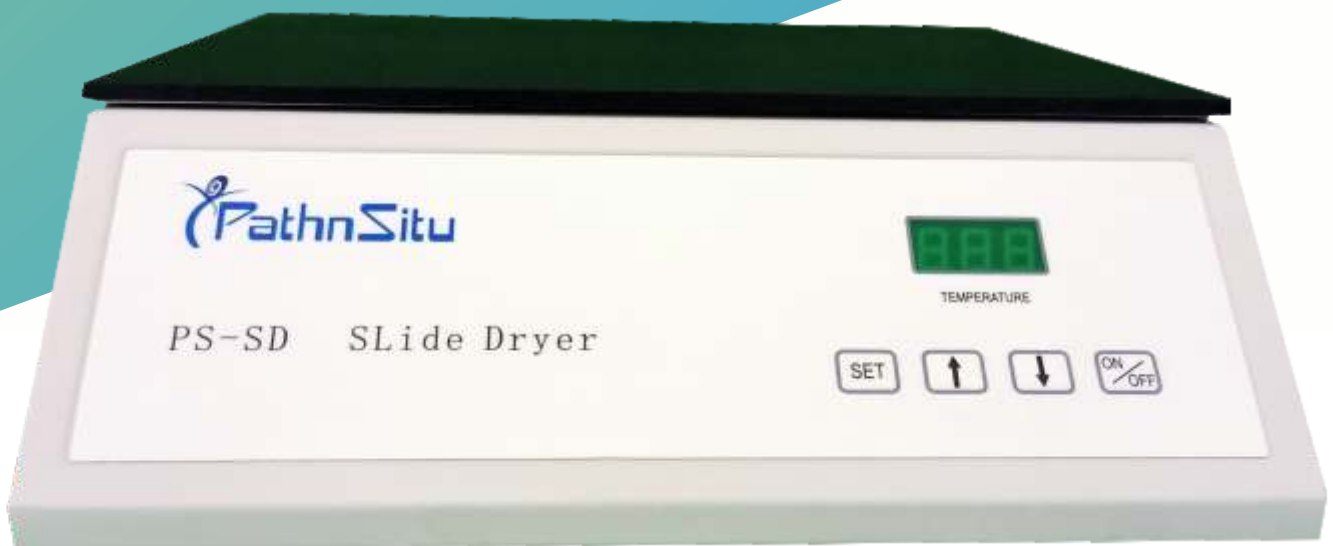


SLIDE DRYER

PS-SD



Features

- Provides quick heat transmission on the plate surface
- Maintains precise and uniform temperature
- User friendly set temperature with LED display
- Black surface of the plate provides clear visibility of the tissue sections
- Wide surface of hotplate can accommodate around 35 slides

Description

Slide Dryer or Slide warming table is generally used for slide baking. The step of slide baking is where the FFPE tissue slides after the sectioning are baked/heated using dry heat to improve the adhesion of tissue sections to the slide.

Technical Specifications

Size of heating plate (in mm)	340 x 260
Temperature set range	Room temperature to 90°C ± 2°C
Voltage	AC 220±10% 50Hz
Power	200W
Dimensions (in mm)	350(L)×340(W)×120(H)
Weight	5kg

Ordering Information

Equipment Name	Catalog No
Slide Dryer/Slide Warming Table	PS-SD

Contact Us

USA

1257, Quarry Ln, Suite150
Pleasanton, CA-94566, USA
Phone: +1 925-218-6939

Email

techsupport@pathnsitu.com
customerservice@pathnsitu.com

www.pathnsitu.com

India

PS Heights, 2nd Floor, Plot# A-28/1/11B/B,
Road#15 Nacharam IDA,
Hyderabad-500 076, Telangana.

Phone: 040-2701-5533 /5544
Fax: 040-2701-5544

